

## MacConkey Agar RS Plate

MP1702

### Intended Use:

Recommended for isolating and differentiating Gram negative enteric bacilli from specimens containing swarming strains of *Proteus* species.

### Composition\*\*

Ingredients	g / L
Peptone	17.000
Proteose peptone	3.000
Lactose	10.000
Bile salts	5.000
Sodium chloride	5.000
Neutral red	0.030
Crystal violet	0.001
Agar	13.500
Final pH ( at 25°C)	7.1±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

### Principle And Interpretation

MacConkey agars are slightly selective and differential plating media mainly used for the detection and isolation of gram-negative organisms from clinical (1,2), dairy(3), food(4,5), water(6), and industrial sources (7). It is also recommended for the selection and recovery of the Enterobacteriaceae and related enteric gram-negative bacilli. These agar media is less selective since the concentration of bile salts, which inhibit gram-positive microorganisms, is low in comparison with other enteric plating media. Other than that this medium is also used for count of coli-aerogenes bacteria in cattle and sheep faeces (1), the count of coli-aerogenes and non-lactose fermenters in poultry carcasses (8), bacterial counts on irradiated canned minced chicken (9).

MacConkey Agar is the earliest selective and differential medium for cultivation of enteric microorganisms from a variety of clinical specimens (10,11). The original medium contains protein, bile salts, sodium chloride and two dyes. The selective action of this medium is attributed to crystal violet and bile salts, which are inhibitory to most species of gram-positive bacteria. Gram-negative bacteria usually grow well on the medium and are differentiated by their ability to ferment lactose. Lactose-fermenting strains grow as red or pink colonies and may be surrounded by a zone of acid precipitated bile. The red colour is due to production of acid from lactose, absorption of neutral red and a subsequent colour change of the dye when the pH of medium falls below 6.8. Lactose non-fermenting strains, such as *Shigella* and *Salmonella* are colourless, transparent and typically do not alter appearance of the medium.

Peptones are sources of nitrogen and other nutrients. Lactose is a fermentable carbohydrate, bile salts and crystal violet are selective agents that inhibit growth of gram-positive organisms. Neutral red is the pH indicator dye.

### Type of specimen

Clinical - faeces, urine and other pathological material, foodstuffs and dairy samples, water samples, pharmaceutical samples.

### Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (1,2).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (3-5).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards(6).

For industrial samples, follow appropriate techniques for sample collection, processing as per guidelines (7).

After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions

In Vitro diagnostic use only. Read the label before opening the container. Wear protective gloves/ protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens

and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

## Limitations

1. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium
2. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
3. It is recommended to store the plates at 24-30°C to avoid minimum condensation.
4. Though the medium is recommended for selective isolation, further biochemical and serological testing must be carried out for further confirmation.
5. The surface of the medium should be dry when inoculated.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Sterile MacConkey Agar RS Plate in 90 mm disposable plates with smooth surface and absence of black particles/cracks/bubbles

### Colour of medium

Orange red coloured medium

### Quantity of medium

25ml of medium in 90mm disposable plate

### pH

6.90-7.30

### Cultural Response

Cultural response was observed after an incubation at 30-35°C for 18-72 hours. Recovery rate is considered as 100% for bacteria growth on Soybean Casein Digest Agar.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony
<i>Corynebacterium diphtheriae</i> type gravis	$\geq 10^3$	inhibited	0%	
<i>Shigella flexneri</i> ATCC 12022 (00126*)	50 -100	fair to good	30 -40 %	colourless
<i>Salmonella</i> Paratyphi A ATCC 9150	50 -100	luxuriant	$\geq 50$ %	colourless
<i>Salmonella</i> Abony NCTC 6017 (00029*)	50 -100	luxuriant	$\geq 50$ %	colourless
## <i>Proteus hauseri</i> ATCC 13315	50 -100	luxuriant	$\geq 50$ %	colourless
<i>Salmonella</i> Typhi ATCC 6539	50 -100	luxuriant	$\geq 50$ %	colourless
<i>Staphylococcus epidermidis</i> ATCC 12228 (00036*)	$\geq 10^3$	inhibited	0%	
<i>Escherichia coli</i> ATCC 8739 (00012*)	50 -100	luxuriant	$\geq 50$ %	pink-red with bile precipitate
<i>Staphylococcus aureus</i> subsp.aureus ATCC 6538 (00032*)	$\geq 10^3$	inhibited	0%	
<i>Salmonella</i> Paratyphi B ATCC 8759	50 -100	luxuriant	$\geq 50$ %	colourless
<i>Escherichia coli</i> ATCC 25922 (00013*)	50 -100	luxuriant	$\geq 50$ %	pink to red with bile precipitate
# <i>Klebsiella aerogenes</i> ATCC 13048 (00175*)	50 -100	luxuriant	$\geq 50$ %	pink to red

<i>Salmonella</i> Typhimurium ATCC 14028 (00031*)	50 -100	luxuriant	>=50 %	colourless
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	50 -100	none - poor	<=10 %	colourless to pale pink
<i>Salmonella</i> Enteritidis ATCC 13076 (00030*)	50 -100	luxuriant	>=50 %	colourless
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	>=10 <sup>3</sup>	inhibited	0%	

Key :- \* Corresponding WDCM numbers

# Formerly known as *Enterobacter aerogenes*

## Formerly known as *Proteus vulgaris*

## Storage and Shelf Life

On receipt store between 20-30°C Use before expiry date on the label. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (1,2).

## Reference

1. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
2. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
3. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.
4. FDA Bacteriological Analytical Manual, 2005, 18th Ed., AOAC, Washington, D.C.
5. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
6. Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.
7. Williams, (Ed.), 2005, Official Methods of Analysis of the Association of Official Analytical Chemists, 19th Ed., AOAC, Washington, D.C
8. Barnes Ella M. and Shrimpton D. H., 1957, J. Appl. Bacteriol., 20(2), 273-285.
9. Thornley Margaret J., 1957, J. Appl. Bacteriol., 20(2), 273-285.
10. MacConkey A., 1905, J. Hyg., 5:333.
11. MacConkey A., 1900, The Lancet, ii:20.


Revision : 01/ 2024



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India



*In vitro* diagnostic  
medical device

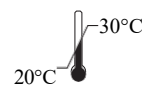

Storage temperature



CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
[www.cepartner4u.eu](http://www.cepartner4u.eu)



CE Marking



Do not use if  
package is damaged

#### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.