

## OFBL Agar Plate (Oxidation Fermentation Polymyxin Bacitracin Lactose Agar Plate)

MP1811

### Intended Use:

With Polymyxin and Bacitracin is recommended for the selective isolation of *Burkholderia cepacia* from clinical specimens as well as non-clinical samples.

### Composition\*\*

<b>Ingredients</b>	<b>g / L</b>
Tryptone	2.000
Dipotassium hydrogen phosphate	0.300
Sodium chloride	5.000
Lactose	10.000
Bromothymol blue	0.030
Agar	15.000
<b>OFBBL Selective Supplement</b>	<b>FD269</b>
Bacitracin	200units
Polymyxin B sulphate	300,000units
Final pH ( at 25°C)	6.8±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

### Principle And Interpretation

*Burkholderia cepacia* is an opportunistic pathogen generally associated with nosocomial infections. Due to its ability to survive for extended period of time in hostile environments, it is found in such widely varied and inhibitory items such as equipment, medications, mouthwash and disinfectants. Nosocomial infections caused by this organism include bacteremia, urinary infections, and respiratory infections. However, the most serious implication is when identified in patients with Cystic Fibrosis (CF) Patients with cystic fibrosis (CF) have a predisposition for infection and infected patients, if untreated, show a rapid decline in lung function, frequent bacteremia, and death due to lung failure. It is also reported to be a primary cause of bacteremia, pneumonia, and death in the Chronic granulomatous disease CGD patient population (1,2). Therefore, it is critical that isolation and proper identification be fast and accurate. *Burkholderia cepacia* Agar Plate (MP1640) as well as OFBL Agar Plate (MP1811) is recommended for isolation of *Burkholderia cepacia* from clinical specimens (1,3,4).

OFBBL medium is a modified version of OF (Oxidation Fermentation) basal medium developed by a group of researchers, Welch et al. in 1987 (5). They demonstrated that OFBBL agar resulted in improved recovery and isolation of *Burkholderia cepacia* when compared to other selective mediums such as MacConkey agar, XLD agar and various blood blood agars. OFBBL Agar contains the Tryptone which provides necessary nitrogenous compounds and lactose serves as carbohydrate source. Lactose is readily utilized by *Burkholderia cepacia*. The fermentation of lactose results in the release of acid end-products which is detected by the pH indicator, bromothymol blue, present in the medium. When sufficient acid is produced the medium changes from green to yellow providing the colonies their yellow coloration. Dipotassium hydrogen phosphate in the medium buffers the medium well. Sodium chloride helps to maintain osmotic balance. The selectivity of the medium owes itself to the presence of the antibiotics polymyxin B and bacitracin together these antibiotics provide good suppression of the bacterial flora present in respiratory secretions and sputum (for the inhibition of gram-positive organisms and *Neisseria*) (5).

Appropriate patient samples for testing include sputum, bronchial washings, and pharyngeal swabs. Using a direct or diluted inoculum from the sample, a four-quadrant streak may be performed to obtain well-isolated colonies. If the sample is collected on a swab, the swab may be rolled over a small area near the edge of the plate and further streaked for isolation using a sterile loop. It is recommended to incubate at 30-35°C and examine the plates daily upto five days. before discarding. Generally colonies develop within 48-72 hours.

Typically, *Burkholderia cepacia* colonies appear as yellow colonies with yellow halos. Most strains will grow in 48 hours but some strains may require up to 5 days for the color development. Most other microorganisms are inhibited on this medium. Other similar bacteria such as *Pseudomonas aeruginosa* may also grow as yellow colonies and further testing be performed to differentiate these bacteria from *Burkholderia cepacia*. *Burkholderia gladioli* which has been shown to occur in respiratory tract specimens of CF patients will grow on OFPBL Agar and may resemble *B.cepacia* (6). Any growth is considered a positive result and additional biochemical and/or serological tests should be performed on isolated colonies from pure culture (1). The yellowing of the medium signifies carbohydrate fermentation and some rare, atypical strains of *Burkholderia cepacia* may not produce this reaction.

## Type of specimen

Clinical samples - bronchial washings, and pharyngeal swabs

## Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (7,8). After use, contaminated materials must be sterilized by autoclaving before discarding.

## Warning and Precautions

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

## Limitations

1. Similar bacteria such as *Pseudomonas aeruginosa* may also grow as yellow colonies and further testing must be performed to differentiate these bacteria from *Burkholderia cepacia*.
2. Further biochemical and serological tests must be carried out for complete identification.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Sterile O.F.P.B.L Agar in 90 mm disposable plates with smooth surface and absence of black particles/ cracks/ bubbles.

### Colour of medium

Green coloured medium.

### Quantity of medium

25ml of medium in 90mm of disposable plates

### pH

6.60-7.00

### Sterility Check

Passes release criteria

### Cultural Response

Cultural characteristics observed after an incubation at 30-35°C for 48 - 72 hours. Examine upto five days.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony
<i>Burkholderia cepacia</i> ATCC 25416	50-100	good-luxuriant	>=50%	yellow w/ yellow halo
<i>Escherichia coli</i> ATCC 25922 (00013*)	>=10 <sup>3</sup>	inhibited	0%	
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	>=10 <sup>3</sup>	inhibited	0%	

Key : (\*) Corresponding WDCM numbers.

## Storage and Shelf Life

On receipt store between 2-8°C. Use before expiry date on the label. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (7,8).

## Reference

1. Gilligan, P. H. and P. Vandamme, 2003. Misc. Gram Negative Bacteria, pp 729-748. In Murray, P. R., et al., Manual of Clinical Microbiology, 8th ed., American Society for Microbiology, Washington D.C., 2003.
2. O'Neil, K.M., et al. 1986. J. Pediatr. 108: 940-942.
3. Carson, L.A. et al. 1988. J. Clin. Microbiol. 25:1730-1734.
4. Tablan, O.C., et al. 1987. J. Clin. Microbiol. 25:485-487.
5. Welch DF, Muszynski M. J, Pai CH, Marcon MJ, Hribar MM, Gilligan PH, Matsen JM, Ahlin PA, Hilman BC, Chartrand SA. 1987. J. Clin. Microbiol; 25:1730-4.
6. Christenson, J.C., et al. 1989. J. Clin. Microbiol. 27: 270-273..
7. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
8. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision: 02/2024



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India



**IVD** *In vitro diagnostic  
medical device*

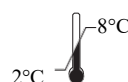

**Storage temperature**



CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
www.cepartner4u.eu



**CE Marking**



**Do not use if  
package is damaged**

### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.