

## HiCrome™ Coliform Agar Plate, Modified

MP1832

### Intended Use:

Recommended for the simultaneous detection of *Escherichia coli* and total coliforms in water and food samples.

### Composition\*\*

Ingredients	Gms / Litre
Peptone, special	8.000
Sodium chloride	1.000
Yeast extract	3.000
Potassium dihydrogen phosphate	0.200
Dipotassium hydrogen phosphate	0.600
Bile Salts	0.800
Magnesium sulphate	0.200
Chromogenic mixture	0.200
Agar	10.000
Final pH ( at 25°C)	7.2±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

### Principle And Interpretation

HiCrome™ Coliform Agar Modified is a selective medium recommended for the simultaneous detection of *E.coli* and thermo-tolerant coliforms in water and food samples (1). Peptone special and yeast extract provide carbon, nitrogen compounds, long chain amino acids, vitamins and essential growth nutrients to the organisms. The phosphates buffer the medium well. Magnesium sulphate helps colour development. Bile salts inhibits gram-positive organisms. Sodium chloride maintains osmotic balance. The chromogenic mixture contains two chromogenic substrates, which enables the detection of two specific enzymes,  $\beta$ -galactosidase and  $\beta$ -glucuronidase.  $\beta$ -galactosidase produced by coliforms cleaves one chromogen, resulting in the pink colouration of coliform colonies. The enzyme  $\beta$ -glucuronidase produced by *E.coli*, cleaves X glucuronide. *E.coli* forms dark blue to violet coloured colonies due to cleavage of both the chromogens (2,3,4). *E.coli* strains that are  $\beta$ -glucuronidase negative (serotype O157:H7) produce pink coloured colonies. Other gram negative bacteria able to grow at (44±0.5 °C) produce white or colourless colonies. Transfer 1 ml of product to analyze and its tenfold dilutions to sterile Petri plates. Pour 12 ml of medium, mix well and allow to solidify. Overlay with 4 ml of medium, allow to solidify and incubate at 43-45°C for 18-24 hours.

### Type of specimen

Clinical sample-urine, faeces; Food and Dairy samples; Water samples.

### Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (5,6).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (7-9).

For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (10). After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions

In Vitro diagnostic use. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

### Limitations

1.  $\beta$ -glucuronidase is present in 97% of *E.coli* strains, however few *E.coli* may be negative.
2. Individual organisms differ in their growth requirement and may show variable growth patterns on the medium.
3. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.

Please refer disclaimer Overleaf.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

### Quality Control

#### Appearance

Sterile HiCrome™ Coliform Agar, Modified in 90mm disposable plates with smooth surface & absence of black particles/ cracks/bubbles.

#### Colour of medium

Light yellow coloured medium

#### Quantity of medium

25 ml of medium in 90 mm disposable plates.

#### pH

7.00-7.40

#### Sterility check

Passes release criteria

#### Cultural Response

Cultural characteristics observed after incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of Colony
<i>Salmonella</i> Enteritidis ATCC 13076 (00030*)	50-100	good	40-50%	colourless
<i>Enterobacter cloacae</i> ATCC 23355	50-100	good-luxuriant	≥50%	pink
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	≥10 <sup>4</sup>	inhibited	0%	
<i>Klebsiella pneumoniae</i> ATCC 13883 (00097*)	50-100	good-luxuriant	≥50%	light pink
<i>Staphylococcus aureus</i> ATCC 25923 (00034*)	≥10 <sup>4</sup>	inhibited	0%	
<i>Escherichia coli</i> ATCC 10536	50-100	good-luxuriant	≥50%	dark blue/violet
<i>Escherichia coli</i> ATCC 25922 (00013*)	50-100	good-luxuriant	≥50%	dark blue/violet

Key : \*Corresponding WDCM numbers.

### Storage and Shelf Life

On receipt store between 2-8°C. Use before expiry date on the label. Product performance is best if used within stated expiry period.

### Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (5,6).

### Reference

1. Manafi M. and Kneifel W., 1989, Zentralbl. Hyg., 189:225.
2. Frampton E. W., Restaino L. and Blaszkowski N., 1988, J. Food Prot., 51:402.
3. Kilian M. and Bülow P., 1976, Acta. Pathol. Microbiol. Scand., Sect. B, 84:245.
4. LeMinor L. and Hamida F., 1962, Ann. Inst. Pasteur (Paris), 102:267.
5. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
6. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
7. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.

8. Salfinger Y., and Tortorello M.L. 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
9. Wehr H.M. and Frank J.H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.
10. Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC: APHA Press; 2023.

Revision : 01/2023



HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India




CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
www.cepartner4u.eu



**In vitro diagnostic  
medical device**



**CE Marking**

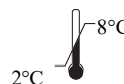

**Storage temperature**



**Do not use if  
package is damaged**

#### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.