

Sabouraud Dextrose Agar Plate w/Gentamicin

MP5386

Intended Use:

Recommended for the selective isolation of fungi from clinical and non-clinical specimens.

Composition**

Ingredients	g / L
Dextrose (Glucose)	40.000
Mycological, peptone	10.000
Gentamicin	0.0400
Agar	15.000
Final pH (at 25°C)	5.6±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Either streak, inoculate or surface spread the test inoculum (50-100 CFU) aseptically on the plate.

Principle And Interpretation

Sabouraud Dextrose Agar is Carliers modification (1) of the formulation described by Sabouraud (2) for the cultivation of fungi (yeasts, moulds), particularly useful for the fungi associated with skin infections. This medium is also employed to determine microbial contamination in food, cosmetics, and clinical specimens (3).

Mycological peptone provides nitrogenous compounds. Dextrose (Glucose) provides an energy source. High dextrose concentration and low pH favours fungal growth and inhibits contaminating bacteria from test samples (4). Gentamicin inhibits most of the bacterial contaminants. Some pathogenic fungi may produce infective spores which are easily dispersed in air, so examination should be carried out in safety cabinet.

Type of specimen

Clinical samples - Skin scraping, nail scraping, etc.; Food samples

Specimen Collection and Handling:

For food samples, follow appropriate techniques for sample collection and processing as per guidelines (5). For clinical samples follow appropriate techniques for handling specimens as per established guidelines (6,7). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Individual strain of a microorganism may have unique growth requirements with respect to nutrients and physical conditions. Based on which the growth pattern of each varies on a medium and some even may display significant delay.
2. Some pathogenic fungi may produce infective spores which are easily dispersed in air, so examination should be carried out in safety cabinet.
3. Each lot of the medium has been tested for the organisms specified on the COA. It is recommended to users to validate the medium for any specific microorganism other than mentioned in the COA based on the user's unique requirement.
4. Further biochemical tests should be carried out for complete identification.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Sterile Sabouraud Dextrose Agar w/ Gentamicin in 90mm disposable plate with smooth surface and absence of black particles/cracks/bubbles.

Colour of Medium

Light amber coloured medium.

Quantity of Medium

25ml of medium in 90mm disposable plates

pH

5.40-5.80

Sterility Check

Passes release criteria

Cultural Response

Cultural characteristics observed after incubation at 22-28 °C for 48-72 hours.

Organism	Inoculum (CFU)	Growth	Recovery
<i>Candida albicans</i> ATCC 10231 (00054*)	50 -100	good to luxuriant (white colonies)	>=50 %
# <i>Aspergillus brasiliensis</i> ATCC 16404 (00053*)	50 -100	good to luxuriant	>=50 %
<i>Saccharomyces cerevisiae</i> ATCC 9763 (00058*)	50 -100	good to luxuriant	>=50 %
<i>Trichophyton interdigitale</i> ATCC 9533	50 -100	good to luxuriant	>=50 %
<i>Escherichia coli</i> ATCC 25922 (00013*)	50 -100	none-poor	0 -10 %
<i>Escherichia coli</i> ATCC 8739 (00012*)	50 -100	none-poor	0 -10 %
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923 (00034*)	50 -100	none-poor	0 -10 %
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 6538 (00032*)	50 -100	none-poor	0 -10 %
<i>Pseudomonas aeruginosa</i> ATCC 27853 (00025*)	50 -100	none-poor	0 -10 %
^ <i>Pseudomonas paraaeruginosa</i> ATCC 9027 (00026*)	50 -100	none-poor	0 -10 %

Key : (*) Corresponding WDCM numbers, (*)-Formerly known as *Aspergillus niger*,
(^) Formerly known as *Pseudomonas aeruginosa* , (\$) Formerly known as *Trichophyton mentagrophytes*.

Storage and Shelf Life

Store between 2-8°C. Use before expiry date on the label. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (6,7).

Reference

1. Carlier G. I. M., 1948, Brit. J. Derm. Syph., 60:61.
2. Sabouraud K., 1892, Ann. Dermatol. Syphilol, 3:1061.
3. Bacteriological Analytical Manual, 8th Edition, Revision A, 1998. AOAC, Washington D.C.
4. Murray PR, Baren EJ, Jorgensen JH, Pfaller MA, Tenover FC, Tenover FC (editors) 2003, Manual of clinical Microbiology, 8th ed., ASM, Washington, D.C.
5. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
6. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
7. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W. (2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.


Revision : 02/2024



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India



IVD *In vitro* diagnostic
medical device

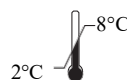

Storage temperature



CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



CE Marking



**Do not use if
package is damaged**

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.