

## TCBS HiVeg<sup>®</sup> Agar

MV189

### Intended Use

Recommended for the selective isolation and cultivation of *Vibrio cholerae* and other enteropathogenic *Vibrio's* causing food poisoning from clinical and food specimen.

### Composition

Ingredients	g/ L
HiVeg <sup>®</sup> peptone No.3	15.000
Yeast extract	6.000
Sodium thiosulphate	10.000
Sodium citrate	10.000
Synthetic detergent No. II	2.000
Sucrose	20.000
Sodium chloride	10.000
Ferric citrate	1.000
Bromo thymol blue	0.040
Thymol blue	0.040
Agar	15.000
Final pH ( at 25°C)	8.6±0.2

\*\*Formula adjusted, standardized to suit performance parameters

### Directions

Suspend 89.08 grams in 1000 ml purified/ distilled water. Heat to boiling to dissolve the medium completely. DO NOT AUTOCLAVE. Cool to 45-50°C. Mix well and pour into sterile Petri plates.

### Principle And Interpretation

TCBS Agar was developed by Kobayashi et al (1), who modified the selective medium of Nakanishi (2). Although this medium was originally designed for the isolation of *V.cholerae* and *V. parahaemolyticus*, most *Vibrio's* grow to healthy large colonies with many different colonial morphologies. TCBS Agar is also recommended by APHA for the selective isolation of *V. cholerae* and *V. parahaemolyticus* (3,4). Enrichment in Alkaline HiVeg<sup>®</sup> Peptone Water (MV618), followed by isolation on TCBS HiVeg<sup>®</sup> Agar is routinely used for isolation of *V.cholerae* (5,6,7). TCBS HiVeg<sup>®</sup> Agar is prepared by using vegetable peptones in place of animal based peptones which make the media free of BSE/TSE risks. HiVeg<sup>®</sup> peptone No.3 and yeast extract provide nitrogenous compounds, vitamin B complex and other essential growth nutrients. Synthetic detergent No. II and sodium citrate inhibit gram-positive bacteria and coliforms (8). Sodium thiosulphate serves as a good source of sulphur, which in combination with ferric citrate detects the production of hydrogen sulphide. For the metabolism of *Vibrio's*, sucrose is added as a fermentable carbohydrate. *Vibrio* that is able to utilize sucrose will form yellow colonies. Bromothymol blue and thymol blue are the pH indicators. The alkaline pH of the medium improves the recovery of *V.cholerae*. Strains of *V. cholerae* produce yellow colonies on TCBS HiVeg<sup>®</sup> Agar because of fermentation of sucrose. *V. alginolyticus* also produce yellow colonies. *V.parahaemolyticus* is a sucrose non-fermenting organism and therefore produces blue-green colonies, as does *V.vulnificus*.

### Type of specimen

Clinical : faeces, etc; Food samples; Water samples.

### Specimen Collection and Handling

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (7,8).

For food and dairy samples, follow appropriate techniques for sample collection and processing as per guidelines (4,9,10). For water samples, follow appropriate techniques for sample collection, processing as per guidelines and local standards (3). After use, contaminated materials must be sterilized by autoclaving before discarding.

### Warning and Precautions

In Vitro diagnostic use. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

## Limitations

1. The medium should be inoculated heavily with faecal specimens because growth of few species may be inhibited on the medium due to fermentation of sucrose and accumulation of acids.
2. However, occasional isolates of *Pseudomonas* and *Aeromonas* may also form blue green colonies on TCBS HiVeg® Agar (11).
3. *Proteus* species that are sucrose-fermenters may form yellow colonies (11).
4. TCBS HiVeg® Agar is not a suitable medium for oxidase testing of *Vibrio* species (12).
5. A few strains of *V. cholerae* may appear green or colourless on TCBS Agar due to delayed sucrose fermentation (11).
6. TCBS HiVeg® Agar is highly selective for *Vibrio* species. Any H<sub>2</sub>S negative colony of TCBS HiVeg® Agar can be considered presumptive positive for *Vibrio*.
7. Further biochemical and serological tests must be carried out for complete identification.

## Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

## Quality Control

### Appearance

Light yellow to light tan homogeneous free flowing powder

### Gelling

Firm, comparable with 1.5% Agar gel

### Colour and Clarity of prepared medium

Bluish green coloured clear to slightly opalescent gel forms in Petri plates.

### Reaction

Reaction of 8.9% w/v aqueous solution at 25°C. pH : 8.6±0.2

### pH

8.40-8.80

### Cultural Response

Cultural characteristics observed after an incubation at 35-37°C for 18-24 hours.

Organism	Inoculum (CFU)	Growth	Recovery	Colour of colony
<b>Productivity</b>				
<i>Vibrio parahaemolyticus</i> NCTC 10885 (00185*)	50-100	good-luxuriant	≥50%	blue
<i>Vibrio furnissii</i> NCTC 11218 ((00186*)*)	50-100	good-luxuriant	≥50%	greenish yellow
<b>Selectivity</b>				
<i>Escherichia coli</i> ATCC 25922 (00013*)	≥10 <sup>4</sup>	inhibited	0%	
<i>Escherichia coli</i> ATCC 8739(00012*)	≥10 <sup>4</sup>	inhibited	0%	
<i>Escherichia coli</i> ATCC 11775 (00090*)	≥10 <sup>4</sup>	inhibited	0%	
<b>Additional Microbiological Testing</b>				
<i>Vibrio parahaemolyticus</i> ATCC 17802 (00037*)	50-100	good-luxuriant	≥50%	blue
<i>Vibrio cholerae</i> ATCC 15748	50-100	good-luxuriant	≥50%	yellow
<i>Vibrio fluvialis</i> ATCC 33809 (00137*)	50-100	good-luxuriant	≥50%	yellow
<i>Vibrio vulnificus</i> ATCC 29306	50-100	fair - good	≥20%	greenish yellow

## <i>Proteus hauseri</i> ATCC 13315	>=10 <sup>4</sup>	inhibited	0%
<i>Shigella flexneri</i> ATCC 12022 (00126*)	>=10 <sup>4</sup>	inhibited	0%
<i>Enterococcus faecalis</i> ATCC 29212 (00087*)	>=10 <sup>4</sup>	inhibited	0%

Key: (\*) Corresponding WDCM numbers

## Formerly known as *Proteus vulgaris*

## Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 20-30°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

## Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (7,8).

## Reference

1. Kobayashi T., Enomoto S., Sakazaki R., and Kuwahara S., 1963, Jap. J. Bacteriol., 18: 387.
2. Nakanishi Y., 1963, Modern Media 9: 246.
3. Lipps WC, Braun-Howland EB, Baxter TE, eds. Standard methods for the Examination of Water and Wastewater, 24th ed. Washington DC:APHA Press; 2023.
4. Salfinger Y., and Tortorello M.L. Fifth (Ed.), 2015, Compendium of Methods for the Microbiological Examination of Foods, 5th Ed., American Public Health Association, Washington, D.C.
5. Forbes B. A., Sahn A. S. and Weissfeld D. F., 1998, Bailey & Scotts Diagnostic Microbiology, 10th Ed., Mosby, Inc. St. Louis, Mo.
6. Furniss A. L., Lee J. V. and Donovan T. J., 1978, The Vibrios, Public Health Laboratory Service Monograph Series No. 11, Maidstone Public Health Laboratory, H.M.S.O., London, England.
7. Murray P. R., Baron J. H., Pfaller M. A., Jorgensen J. H. and Tenover F. C., (Eds.), 2003, Manual of Clinical Microbiology, 8th Ed., American Society for Microbiology, Washington, D.C.
8. Isenberg, H.D. Clinical Microbiology Procedures Handbook. 2nd Edition.
9. Wehr H. M. and Frank J. H., 2004, Standard Methods for the Microbiological Examination of Dairy Products, 17th Ed., APHA Inc., Washington, D.C.
10. American Public Health Association, Standard Methods for the Examination of Dairy Products, 1978, 14th Ed., Washington D.C.
11. MacFaddin J. F., 1985, Media for Isolation-Cultivation-Identification-Maintenance of Medical Bacteria, Vol. 1, Williams & Wilkins, Baltimore, Md.
12. Morris G. K., Merson M. H., Huq A. K., Kibria A. K. and Black R., 1979, J. Clin. Microbiol., 9:79


Revision : 04/2025




HiMedia Laboratories Pvt. Limited,  
Plot No.C-40, Road No.21Y,  
MIDC, Wagle Industrial Area,  
Thane (W) -400604, MS, India



**In vitro diagnostic  
medical device**

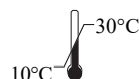

**Storage temperature**

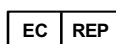

CEpartner4U, Esdoornlaan 13,  
3951DB Maarn, NL  
www.cepartner4u.eu



**CE Marking**



**Do not use if  
package is damaged**

### Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.