

Listeria Selective HiVeg™ Broth Base

MV889

Intended Use:

Recommended for selective isolation and cultivation of *Listeria monocytogenes* from clinical specimens.

Composition**

Ingredients	Gms / Litre
HiVeg™ hydrolysate	17.000
Soya peptone	3.000
Yeast extract	6.000
Sodium chloride	5.000
Dipotassium hydrogen phosphate	2.500
Dextrose (Glucose)	2.500
Final pH (at 25°C)	7.3±0.2

**Formula adjusted, standardized to suit performance parameters

Directions

Suspend 36.0 grams in 1000 ml purified / distilled water. Heat if necessary to dissolve the medium completely. Sterilize by autoclaving at 15 lbs pressure (121°C) for 15 minutes. Cool to room temperature and aseptically add rehydrated contents of one vial of Listeria Selective Supplement II (FD063) or two vials of Listeria Selective Supplement II (FD063I) as desired. Mix well before dispensing into sterile tubes or flasks as desired.

Principle And Interpretation

Listeria monocytogenes is a short, gram-positive, non spore-forming rod shaped bacterium that appears coccoidal in older cultures. *Listeria* multiplies over a wide range of temperatures from 3°C to 45°C with optimum temperature range of 30°C to 37°C. *L. monocytogenes* has been isolated from numerous environmental sources such as silage, soil, decaying vegetation, sewage, damp earth, straw and faeces (2,9). Detection of *L.monocytogenes* in foods is not difficult. Low numbers of organisms are commonly isolated from raw milk, meat, vegetables, seafood and the food-processing environment. Enrichment procedures are used to isolate low numbers of *L.monocytogenes*. Injured *L. monocytogenes* are sublethally stressed as a result of exposure to heat, freezing or acidic conditions. Sublethally stressed *L. monocytogenes* require resuscitation in a non-selective medium at a temperature favouring repair of the sublethal injury.

Listeria Selective Broth is formulated as per Lovett et al (7) for the selective enrichment of *Listeria* species from milk and milk products and other foods. Listeria Selective Broth is recommended by ISO Committee (3) with a slight modification in the supplement (FD063I). Listeria Selective HiVeg™ Broth Base is prepared by completely replacing animal based peptone with vegetable peptones to avoid BSE/TSE risks associate with animal peptones.

HiVeg™ hydrolysate, soya peptone and yeast extract provide carbon and nitrogen compounds, long chain amino acids, vitamins essential for bacterial metabolism. Dextrose is the energy source. The medium is rendered selective by addition of selective supplement. Cycloheximide inhibits the growth of saprophytic fungi. Nalidixic acid inhibits growth of gram-negative organisms and acriflavin suppresses gram-positive microorganisms (6,8). Acriflavin is an acridinic derivative with bacteriostatic properties towards many gram-positive bacteria and a weak fungicidal activity.

For enrichment, 25 grams or 25 ml sample is added to 225 ml medium in a stomacher bag. Homogenize the material if required. Incubation is carried out at 30°C for upto 7 days. Ajello et al (1) showed that incubation period of 7 days allows better recovery of environmentally stressed *Listeria* from milk and milk products. The enrichment broth is further subcultured on Listeria Selective HiVeg™ Agar (MV567) after 1, 2 and 7 days.

Listeria monocytogenes is a highly pathogenic organism and therefore proper precautions should be taken while handling them.

Type of specimen

Clinical samples - Body tissues or body fluids

Specimen Collection and Handling:

For clinical samples follow appropriate techniques for handling specimens as per established guidelines (4,5). After use, contaminated materials must be sterilized by autoclaving before discarding.

Warning and Precautions :

In Vitro diagnostic Use only. For professional use only. Read the label before opening the container. Wear protective gloves/protective clothing/eye protection/ face protection. Follow good microbiological lab practices while handling specimens and culture. Standard precautions as per established guidelines should be followed while handling clinical specimens. Safety guidelines may be referred in individual safety data sheets.

Limitations :

1. Further biochemical tests are needed for a final identification of the isolated organisms.
2. *Listeria monocytogenes* is a highly pathogenic organism and therefore proper precautions should be taken while handling them.

Performance and Evaluation

Performance of the medium is expected when used as per the direction on the label within the expiry period when stored at recommended temperature.

Quality Control

Appearance

Cream to yellow homogeneous free flowing powder

Colour and Clarity of prepared medium

Fluorescent yellow coloured, clear solution in tubes

Reaction

Reaction of 3.6% w/v aqueous solution at 25°C. pH : 7.3±0.2

pH

7.10-7.50

Cultural Response

Cultural characteristics observed with added Listeria Selective Supplement II (FD063 / FD063I) after an incubation at 30-35°C for 24-48 hours.

Organism	Inoculum (CFU)	Growth
<i>Candida albicans</i> ATCC 10231 (00054*)	≥10 ⁴	inhibited
<i>Escherichia coli</i> ATCC 25922 (00013*)	≥10 ⁴	inhibited
<i>Listeria monocytogenes</i> ATCC 19111 (00020*)	50-100	luxuriant
<i>Listeria monocytogenes</i> ATCC 19112	50-100	luxuriant
<i>Listeria monocytogenes</i> ATCC 19118	50-100	luxuriant
<i>Staphylococcus aureus</i> subsp. <i>aureus</i> ATCC 25923(00034*)	50-100	none-poor

Key : *Corresponding WDCM numbers.

Storage and Shelf Life

Store between 10-30°C in a tightly closed container and the prepared medium at 2-8°C. Use before expiry date on the label. On opening, product should be properly stored dry, after tightly capping the bottle in order to prevent lump formation due to the hygroscopic nature of the product. Improper storage of the product may lead to lump formation. Store in dry ventilated area protected from extremes of temperature and sources of ignition. Seal the container tightly after use. Product performance is best if used within stated expiry period.

Disposal

User must ensure safe disposal by autoclaving and/or incineration of used or unusable preparations of this product. Follow established laboratory procedures in disposing of infectious materials and material that comes into contact with clinical sample must be decontaminated and disposed of in accordance with current laboratory techniques (4,5).

Reference

1. Ajello G., Hayes P. and Fuley J., 1986, Abstracts of the Annual Meeting, ASM, Washington, D.C.
2. Gray M. L., 1960, Science, 132:1767.
3. International Organization for Standardization (ISO), 1993, 10560 Ind. Technical, Corrigendum Cor. 1:1994.

4. Isenberg, H.D. Clinical Microbiology Procedures Handbook 2nd Edition.
5. Jorgensen, J.H., Pfaller, M.A., Carroll, K.C., Funke, G., Landry, M.L., Richter, S.S and Warnock., D.W.(2015) Manual of Clinical Microbiology, 11th Edition. Vol. 1.
6. Lee W. K. and McClain D., 1986, Appl. Environ, Microbiol., 52:1215.
7. Lovette J., Francis D. W. and Hunt J. M., 1987, J. Food Prot., 50:188.
8. McClain D. and Lee W. H., 1988, J. Assoc. off. Anal. Chem., 71:660.
9. Weis J., and Seeliger H. P. R., 1975, Appl. Microbiol. 30:29.

Revision :02/2022



HiMedia Laboratories Pvt. Limited,
Plot No.C-40, Road No.21Y,
MIDC, Wagle Industrial Area,
Thane (W) -400604, MS, India




CEpartner4U, Esdoornlaan 13,
3951DB Maarn, NL
www.cepartner4u.eu



In vitro diagnostic
medical device



CE Marking

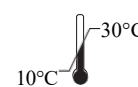

Storage temperature



Do not use if
package is damaged

Disclaimer :

User must ensure suitability of the product(s) in their application prior to use. Products conform solely to the information contained in this and other related HiMedia™ publications. The information contained in this publication is based on our research and development work and is to the best of our knowledge true and accurate. HiMedia™ Laboratories Pvt Ltd reserves the right to make changes to specifications and information related to the products at any time. Products are not intended for human or animal or therapeutic use but for laboratory, diagnostic, research or further manufacturing use only, unless otherwise specified. Statements contained herein should not be considered as a warranty of any kind, expressed or implied, and no liability is accepted for infringement of any patents.